



Innovation & Safety
In Healthcare



History

We had a humble beginning as a company in the year 2007 caring for human lives. We at Lars Medicare, aim our efforts in establishing and consolidating ourselves as a leading manufacturer of medical disposables.

We continue developing special know-how in the design, development, and manufacturing of our product range. The journey to excellence is continuous and we intend to carry forward the spirit and meet the requirements of the medical profession.

Our production facility has all the operations under one Roof



An in-house computerized EO Gas Sterilizer, validated and complying with ISO 11135: 2008 requirements. The efficiency and compliance of every sterilization batch are confirmed and ensured by our in-house government-approved analytical laboratory.

The infrastructure comprises Class 10,000 and Class 100,000 clean rooms with state-of-the-art HVAC systems, HEPA filtration, positive pressure, controlled temperature, and regulated humidity, which provide a controlled environment for manufacturing.

Our modern all-electric and proportional hydraulic molding machines enable us to achieve precision injection molding of components.

Our growth approach & aim is to accomplish a large customer network, both domestic & international, as a result of our directed approach towards customer satisfaction through Quality Product, Competitive Cost & On-Schedule Deliveries. Our operations and activities are accredited with an ISO GMP certification and the majority of our products are CE marked.

We also have in place state-of-the-art assembly & packaging machines, which ensure the consistency in product quality.

We strive to enhance our market penetration through the continuous incorporation of the latest technology.

Our spirited skilled engineering team excels in the area of design, development of products and processes, using the latest CAD system and state-of-the-art software.

Further, we also provide complete engineering solutions and services to support conceptual project assignments, and also offer contract manufacturing support.

We are currently exporting to 85 countries located across Globe: Europe (65%), Middle East (15%), Africa (10%), Latin America (5%) and Asia (5%).

“We strive to enhance our market penetration by adding latest technology in our team and creating a new benchmark of Quality and Durability which will exceed our customer’s expectation”.

Conventional
Products



Conventional

Provein™



Additional Port Cap
For Medication



Tricut Bevel Needle



6% Taper Luer
Lock Connection



Wings With Suture Hole

Availability

With wings and injection

Available Size: 14, 16, 17, 18, 20, 22, 24, 26 G

Technical specifications

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	33	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20
26	Violet	0.4-0.6	19	13



Provein α

With long and broad needle hub for a perfect grip.
Without wings and injection valve.

Available Size: 14, 16, 17, 18, 20, 22, 24, 26 G



Provein β

With wings and without injection valve.

Available Size: 14, 16, 17, 18, 20, 22, 24, 26 G



Propen

With broad needle hub.

Available Size: 14, 16, 17, 18, 20, 22 & 24 G

Available In

PTFE/PUR

Packaging Details

Product:- I.V Cannula	No. of Pieces /Carton	Carton Dimension (Cms)		
		L	W	H
Provein (14G 16G 17G 18G)	1000	70	32	18
Provein (20G 22G 24G 26G)	1000	64	32	18
Provein α (14G 16G 17G 18G 20G 22G 24G 26G)	1000	54	29	18
Provein β (14G 16G 17G 18G 20G 22G 24G 26G)	1000	70	32	18
Pro-Pen (14G 16G 17G 18G 20G 22G 24G)	500	32	29	10

Conventional

3 Way Stop Cock Proway



Advanced Features

- A Three Way Stop Cock facilitates multiple lines through single I.V. access.
- It is made up of a transparent medical grade polycarbonate.
- It supports smooth tap rotation at 360 degrees.
- Can be used with intra-arterial and intra-venous systems.
- A red tap indicates if used for arterial purposes.
- Blue tap can be used for Intra Venous Systems.
- The product is also available in Lipid Resistant.



6% Luer Lock Connection



360° Rotating Tap



With Pegs

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH in mm)		
			L	W	H
3 Way stop cock	..	500	54	29	18

Pro Ext Extension Line



Advanced Features

- This product is available in PUR/PVC free material.
- The product is ideal for administering contrast medium/fluids during imaging procedures that require patient to enter a gantry; e.g. CT and MRI Scans.
- It has High quality kink-resistant tubing with smooth inner surfaces for laminar and non – turbulent flow.

Availability

It is available in various tube lengths based on customer requirements.

Packaging Details:

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Ext. Tube	..	200	38	29	42

Technical Specifications

- The company provides the length of tube between 5cms to 250cms.

Proway - EXT Three-way Stop-cock With Extension Line



Advanced Features

- This product is available in PUR/PVC free material.
- It has High quality Kink-Resistant tubing with smooth inner surface for laminar, non turbulent flow.
- Three way stop cock with extension provides separation of administration and insertion sites.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Proway Ext - 3 Way with Ext line	10/20/30 cms	200	38	29	42

USP: Available in customised length.

Conventional

IV Set
Pro-Inject



Bulb Latex



Air Vent With Collar/wings



Smooth Flow Controller

Infusion Gravity Set (I V SET) Pro -Venex



About Product

- The drip Chamber has extra care wings – to support grip and insertion.
- Non Kinkable / Flexible tube
- Drop Option available 20/60 drop/ML
- Vented / Non Vented – Variation Available.

USP : Customized packaging & length versatility

Infusion Set For Oncology Maroon Set



Some Drugs are affected by light of special wave length and that should be protected from direct light.

These products are specially made for Oncology patients, so that these light – sensitive API's and pharmaceutical products should be protected from Light during use/ distribution.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	w	H
I.V Set	..	350	54	32	43

Conventional

Blood Transfusion Set



About Product

Blood components must be filtered during transfusion to remove clots and small clumps of platelets and white blood cells that form during collection and storage.

Features

- Sharp Spike
- Non Kink able / Flexible Tube
- Available in 20 Drops
- Option : Single / Double chamber

Pro -Fusion Delta



Bulb Latex



Vented Spike



Flow Controller



Pro -Fusion

Packaging details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
B.T Set	Vent/Non-Vent	400	53.5	32	42.5

Measured Volume Set For Fluid Volmed Delta Pro



Advanced features

- The Measured Volume Set (M.V Set) is Cylindrical in shape, which has a flexible, clear and calibrated measured volume chamber.
- Its chamber has an in built Injection port and air inlet with bacteria filter.
- There is an Automatic shut off floating valve provided in chamber which helps in preventing the air trapping and also helps to indicate volume of fluid in the chamber.
- There is a 15 micron fluid filter in the drip chamber.
- The Measured Volume Set delivers approximately 60 drops/ml.
- It has a Soft and kink-resistant PVC tube, which is 150 cm long.

Availability

- Choice of Latex flush bulb/"Y" injection port/ Latex free.
- Choice of packing available: Ribbon pouch/ HM pouch.

Specially designed for photosensitive & Oncology Volmed Maroon



Packaging details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	w	H
M.V Set	151/110/100	100	53.5	32	42.5

Specially Designed For Blood Volmed Vital



Technical Specifications

- Chamber in size: 100/110/150ml.
- Needle size: 23 Gauge.

Conventional

Hypodermic Syringes



Advanced features

- They are made up of a high clarity transparent barrel molded by using premium quality medical grade polypropylene.
- It consists of "Thermo Plastic Elastomer" (TPE).
- The needles are manufactured by using imported Cannula with precise tip sharpness and strong bonding strength with the hub.

Packaging Details

Particular	Types of Primary Packing	Inner Carton Dimension LxWxH(in mm)	Quantity (in mm) Inner Carton	Outer Carton (in mm) Dimension	Quantity in no. Outer Carton
2,3 ml	Poly Pack	170x170x160	100x30	850x355x500	3,000
5 ml	Poly Pack	185x170x200	100x20	880x390x415	2,000
10 ml	Poly Pack	205x150x160	50x20	785x430x345	1,000
20 ml	Poly Pack	270x160x150	25x20	800x550x335	500
50 ml	Poly Pack	260x230x160	25x12	540x480x500	300
2 ml	Blister Pack	380x110x81	100x30	680x390x436	3,000
2ml	Blister Pack	380x110x81	100x36	690x390x520	3,600
2 ml	Blister Pack	380x110x81	100x20	690x390x410	2,000
3 ml	Blister Pack	380x110x81	100x30	680x390x436	3,000
5 ml	Blister Pack	380x130x90	100x25	660x390x485	2,500
10 ml	Blister Pack	380x145x130	100x15	750x390x410	1,500
10 ml	Blister Pack	380x145x130	100x10	750x390x300	1,000
20 ml	Blister Pack	395x165x100	50x10	845x410x220	500
Blister Needle	Blister Pack	200x85x43	100x80	820x360x235	8,000

Availability

- Barrel & Plunger: Polypropylene Silicon: Dow Corning.
- Gasket: TPE (Shore Hardness – 55 Shore A, Tensile Strength > 30 Kg/cm²).
- Hub & Cap: Polypropylene.

Technical Specifications

Size	Type	Needle(G)	Specifications
1 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
2 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
3 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
3 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
5 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
5 ml	Prime Luer Slip	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
10 ml	Prime Luer Slip	21Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
10 ml	Prime Luer Slip	21Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
20 ml	Prime Luer Slip	21Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
20 ml	Prime Luer Slip	21Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
2 ml	Eco Luer Type	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
3 ml	Eco Luer Type	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
5 ml	Eco Luer Type	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
10 ml	Eco Luer Type	21Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe
3 ml	Prime Luer Lock Safety Syringe	24Gx1"	Sterile, Non Toxic, Non Pyrogenic, 3 Piece Syringe

IV Flow Regulator



Advanced features

- It has a latex free Y injection site.
- It is made up of medical grade non toxic and non kink-able, soft PVC tubing which is Phthalate free.
- It has a drum regulator adjustment.
- It also consists of a Luer lock connector with air filter.

Availability

- With /Without Needle free connector

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	w	H
Flow Regulator	..	200	38	29	42

Technical Specifications

- Flow rate: from 5 to 250 ml/hr.

Medulla Spinal Needle



- Made of superfine surgical steel.
- Quincke point is designed to spread the dural fibers and helps reduce post-dural puncture headaches.
- Precision-formed side hole permits directional flow of anesthetic and reduces the possibility of the bevel straddling the dura.
- Optimal sized stylet prevents tissue coring and obstruction of needle lumen
- Sterile for single use only.
- Available in Quincke type model only.

Spinal Needle With Quincke Bevel

Gauge	Colour	O.D (in mm)	Effective length (±1 in MM)
15	Blue-Grey	1.8	38, 58, 75, 88, 100, 110, 120
16	White	1.6	
17	Red-Purple	1.4	
18	Pink	1.2	
19	Cream	1.1	
20	Yellow	0.9	
21	Dark Green	0.8	
22	Black	0.7	
23	Deep Blue	0.6	
24	Purple	0.55	
25	Orange	0.5	
26	Brown	0.45	
27	Med-Grey	0.4	

Packing Details

- 50/100 pcs per inner carton.
- Master Carton dimension :46x43x35cm.

Latex Surgical Gloves Hand Care



Advanced features

- It has high elasticity and strength.
- It is made up of high quality natural latex.
- Sterilization is being done by Ethylene Oxide.

Availability

- The option of Powdered or powdered free available as per requirement
- Available sizes: 6, 6.5, 7, 7.5, 8, 8.5.
- Shelf Life: 5 Years from the date of manufacturing.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Latex Surgical Gloves	..	500	68	24	44

Technical Specifications

- Moisture content: below 0.3% per glove.

Physical Dimension (for 18 Gms)

Size 18 Gms	Length (mm) Min.	Palm Width (mm)	Thickness (mm)		
			Cuff (Min)	Palm (Min)	Finger (Min)
5.5	250	70 ± 6	0.11mm	0.14mm	0.16mm
6.0		76 ± 6			
6.5		83 ± 6			
7.0	270	89 ± 6			
7.5		95 ± 6			
8.0		102 ± 6			
8.5		108 ± 6			

Physical Properties

Characteristics	Before Ageing	After Ageing 7±2°c for 168 hrs
ASTM D 3577:2009, IS 13422:1992		
Tensile Strength (Mpa) min.	24	18
Ultimate Elongation (%) min.	750	560
Stress at 500% Elongation (Mpa) Max.	5.5	NA
EN 455-2: 2009		
Minimum force at break	9.0 N	9.0 N

Safety

Products



Safety

Pro-Inject Safe



Advanced features

- The product is made up of medical grade Non Toxic, Non kink-able, soft PVC.
- It also consists of a “Priming Cap”.
- The Drip Set comprises of a Drip Chamber and a Drip Tube.
- The product consists of a Spike with Air vent (Non Vented also available).
- It is Latex free.
- It also comprises of a Flow regulator to monitor flow of I.V. Fluid passing through drip set.

Availability

- It is available with Bulb Latex / Tube Latex.
- The product consists of a Spike with Air vent (Non Vented also available).



Collared Wings



Rotating Male Luer Lock



Y-site /additional Injection Site



Auto Stop Filter

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
I.V.Set Safe	-	400	54	32	43

Technical Specifications

- Our company provides the drip chamber of length 17mm
- The tube length provided: 150cms(but may change as per demand).

I.V Set Air Stop Filter



Advanced features

- The product is made up of medical grade Non Toxic, Non kink-able, soft PVC.
- It consists of a “ Air Stop Filter” feature.
- The Drip Set comprises of a Vented Drip Chamber and sharp spike
- It is Latex free.

Technical specifications

- Our company provides the drip chamber of length 17mm
- The tube length provided: 150cms(Can be customised).

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
I.V Set Air Stop Filter	-	400	54	32	43

I.V Set with Needle -Less Connector



Advanced features

- The product is made up of medical grade Non Toxic, Non kink-able, soft PVC.
- The product consists of a Needle-Less Connector, additional port for medication.
- It is Latex free.
- The Drip Set comprises of a Vented Chamber and sharp spike.
- It also comprises of a Flow regulator to monitor flow of Fluid Controller.
- It safeguard from needle stick injury as the needle-less connector can be used without needle.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
I.V.Set Safe	-	400	54	32	43

Technical Specification

- Our company provides the drip chamber of length 17mm
- The tube length provided: 150cms(Can be customised).

Pro-Flow Delta (with Needle Less Connector)



Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH in mm)		
			L	W	H
Flow Regulator with Needle-less connector	-	200	38	29	42

Technical Specification

Flow rate: from 5 to 250 ml/hr

- *Gives accurate flow.
- *Prevents from needle-stick injuries.

Safety

Provein™ Safe




Port Cap
(with Easy Open
And Close System)



Identification By
Colour Code At Teflon Hub



6% Taper Luer Lock
Connection

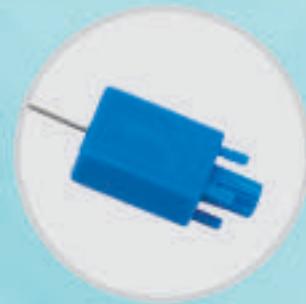

Safety System

Availability

- With wings and injection valve.
- Available Size: 14, 16, 17, 18, 20, 22, 24 G
- Available in PTFE/PUR

Technical Specifications

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	33	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20



Long Projection Body

Provein Safe α

- With long and broad needle hub for a perfect grip.
- Without wings and injection valve.
- Available Size: 14, 16, 17, 18, 20, 22, 24, 26 G




Provein Safe β

- With wings and without injection valve.
- Available Size: 14, 16, 17, 18, 20, 22, 24 G

Packaging details

Product	Size	No. of Pieces/ Carton	Carton Dimensions (LxWxH in mm)		
			L	W	H
Safety I.V Cannula/ Catheter	14,16,17,18,20, 22,24 G	500	45	36	18

Safety

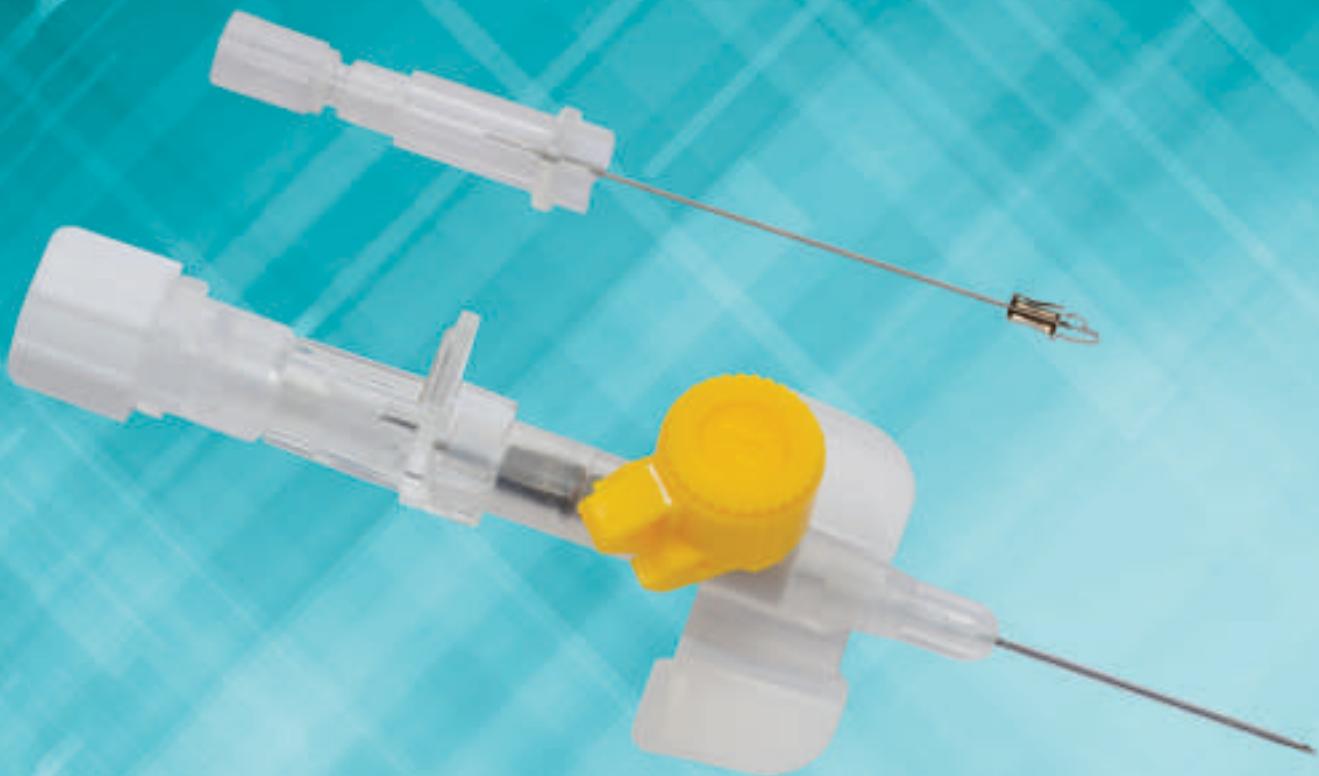

Safe Inject

- With wings and with injection valve.
- Available Size: 14, 16, 17, 18, 20, 22, 24G



Safety Clip



With Wings And Port
Cap For Medication



Safe Cath



Long Projection Body

- * Transparent Hub Cover.
- * Available in PUR/PTFE
- * Sizes : 14, 16, 17, 18, 20, 22, 24G



Flash Back Chamber



Safe Cath Pro

Packaging details

Product	Size	No. of Pieces/ Carton	Carton Dimensions (LxWxH in mm)		
			L	W	H
Safety I.V Cannula/ Catheter	14,16,17,18,20, 22,24 G	500	45	36	18

Technical Specifications

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	33	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

DEHP FREE

Products



Extension Line Pro Ext



Advanced features

- The product is DEHP free.
- The product is ideal for administering contrast medium/fluids during imaging procedures that require patient to enter a gantry; e.g. CT and MRI Scans.
- It has High quality kink-resistant tubing with smooth inner surfaces for laminar and non – turbulent flow.

Technical Specifications

- The company provides the length of tube between 5cms to 250cms.

Availability

- It is available in various tube lengths based on customer requirements.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Ext. Tube	..	200	38	29	42

Three-way Stop-cock With Extension Line Proway Ext



Advanced features

- The product is DEHP free.
- It has High quality Kink-Resistant tubing with smooth inner surface for laminar, non turbulent flow.
- Three way stop cock with extension provides separation of administration and insertion sites.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Proway Ext - 3 Way with Ext line	10/20/30 cms	200	38	29	42

Technical Specifications

- The company provides the length of tube between 5cms to 250cms.

IV SET Pro- Venex



Advanced features

- The product is DEHP free.
- The Drip Set comprises of a Drip Chamber and a Drip Tube.
- It also comprises of a Flow regulator to monitor flow of I.V. Fluid passing through drip set.

Respiratory

Products



HME Filter



Advanced Features

- HME filters are passive humidification systems that retain heat and moisture from the patient’s exhaled gas via a hygroscopic filter.
- The product is designed according to the ISO standards that help in facilitating the connection to different ventilators and anesthesia systems.
- The Condensation surface with hygroscopic coating ensures effective moisture exchange and heat retention thus, minimizing the moisture and heat loss of the patient.

Availability

- It comes with a CO₂ monitoring port.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
HME Filter	..	100	38	29	42

Nebulizer With Mask



Advanced Features

- The set includes a Nebulizer Mask, a Nebulizer chamber and a 200 cms long star lumen tube.
- It has a Nebulization rate of 3 cc in 10 minutes.
- It is a soft mask that ensures the patient’s comfort.
- It has an integrated nose clip and elastic band.
- It has starlumen tubing.
- It has a funnel shaped connector at both ends of tube.

Availability

- Tubing length available: 200 cm.
- The product is available in Adult and Pediatric sizes.
- The product is available in Sterile and Non Sterile.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	w	H
Nebulizer with Mask	Adult/ pediatric	100	54	30	36

Lung Exerciser-3 Ball Prociser



Advanced Features

- The transparent and calibrated chambers consist of color coded triple balls allows visual calibration and estimation of the patient’s lungs fitness progress.
- The product is also advised during post surgery rehabilitation.

Availability:

The product has a sleek and stable design with soft corrugated tube and mouth piece.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH)		
			L	W	H
Lung Exerciser(PROCISER)	..	28	54	30	36

Technical Specifications

- Calibrated and clear three stage chambers for efficient exercise.
- It also has an accommodative groove for tubing when not in use.

Respiratory

Oxygen Mask



Advanced Features

- It has a light weight design.
- It consists of an integrated nose clip and elastic band.
- It consists of star lumen tubing
- It has a funnel shaped connector at both ends.

Availability

- Tubing length available: 200 cm.
- It is available in Adult and Pediatric sizes.
- This product is for a single patient use.
- The product is available in Sterile and Non Sterile.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH in mm)		
			L	W	H
Oxygen Mask	Pediatric/ Adult	200	57	37	38

Technical Specifications

- Integrated nose clip and elastic band for proper mask positioning.
- Star Lumen tubing to ensure continuous flow of oxygen.

Oxy- Cath



Advanced Features

- It consists of twin nasal tips.
- These nasal tips are fixed over soft & kink-resistant loop of PVC tube.
- It is connected with a 200 cms long Star Lumen tube.

Availability

- The product is available in Adult, Pediatric & Neo natal sizes.
- This product is for a single patient use.
- It is available in Sterile and Non Sterile.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimensions (LxWxH in mm)		
			L	w	H
Twin Bore Nasal Catheter	Pediatric/ Adult	200	45	40	48

Technical Specifications

- Twin nasal tips for an equal distribution of oxygen.
- Nasal tips are fixed over soft and kink resistant loop of PVC tube which is adjustable around the head of the patient.
- It has 99.9% filtration efficiency protection against Bacterial/Viral cross infections.
- It supports low dead Space and low flow resistance.

Yankauer Suction Set

Advanced Features

- This product is suitable for pre-operative removal and secretion of body fluids.
- It has soft and kink-resistant ribbed PVC tubing.
- The Yankauer Handle has an open end with 4 lateral eyes to withstand negative pressure and it also ensures continuous and easy suction of fluid.

Availability

- Available in tube length of 210 cm, 250 cm and 300 cm.
- It is available in Sterile and Non Sterile.

Packaging Details

Product	Size	No. of Pieces/ Carton	Carton Dimension (Cms)		
			L	W	H
Yankauer Suction Set	-	125	60	36	42
Yankauer Suction Set	-	150	53	46	44

Technical Specifications

- The Yankauer Handle is available in standard or crown tip, which is either with or without Vacuum thumb control.

Lars Medicare Pvt. Ltd.

Marketing Office : India

📍 136, Avtar Enclave,
2nd Floor, Paschim Vihar,
New Delhi 110063, India
☎ +91 11 45733068
☎ +91 11 45652878
☎ +91 9899703326
✉ mkt@larsmedicare.com
sales@larsmedicare.com

Marketing Office : Canada

📍 555, Industrial Dr, Suite 105,
Milton, ON L9T 5E1, Canada
☎ +1 647 2995115
☎ +1 905 2050305
✉ sameer@larsmedicare.com

Facility at :

📍 Kila No. 16 & 17, Sultanpur,
Opp. Sports Authority of India,
Near Bahalgarh Chowk,
Sonapat - 131021, Haryana, India
☎ +91 8930111001, +91 8930111190
plant@larsmedicare.com
✉ hr@larsmedicare.com

🌐 www.larsmedicare.com

- The description and specification in the brochure are correct at the time of printing.
- In pursuit of continual improvement, we reserve the right to change the design and specification of any of our products without prior notice.
- All photographs, product designs and trademarks represented in this brochure are the sole proprietary intellectual property of Lars Medicare Pvt. Ltd.